



# Atlantis™ WebOrder

Quick start



# Contents

<b>User registration</b> .....	4
<b>Entering an order</b> .....	7
Order page 1 – Overview .....	8
Order page 2 – Items .....	9
Order page 3 – Parallelism .....	11
Order page 4 – Bill & ship .....	12
<b>Order ticket</b> .....	13
<b>Hold for design approval</b> .....	14
<b>Account settings</b> .....	15
<b>Preferences</b> .....	17
<b>Design profiles</b> .....	18
Add preferences .....	19
Enter new profile .....	20
<b>Help</b> .....	21
Emergence width options .....	22
Model type .....	22

User registration

## Welcome page

Enter [www.atlantisweborder.com](http://www.atlantisweborder.com) into your web browser.

Select language from drop down menu.

For registered users, enter your account e-mail address and password and press **Go** button to log in.

New users must register to access Atlantis™ WebOrder.

Click the **Register** link to start the registration process.

**ATLANTIS™**  
WebOrder

**ASTRATECH**  
**DENTAL**  
*get inspired*

**Welcome to Atlantis™ WebOrder!**

If you have an existing account, please enter your e-mail and password and press the "GO" button to begin.  
If you are not a registered user, click on the "Register" link to set up an account

**Login to enter an order**

Email:

Password:

Remember me

**GO**

[Forgot your password?](#)

[Register](#) to use Atlantis™ WebOrder

**Choose your language:**

English

[Customer Service](#) [Astra Tech](#) [Privacy Policy](#) and [Terms of Service](#)

## Register page

You must register to gain access to the Atlantis™ WebOrder site by completing all required information on the **Register** page. This information may be edited at a later time by the account manager on the **My account** page.

**Note:** If **Bill to** or **Ship to** information is different than **Ordering address**, uncheck the box(es) to enter additional information. To complete your registration, the **Terms of service** must be accepted by clicking the check box.

### Register Atlantis™ WebOrder

User

First name:

Last name:

Email:

Ordering Address

Business name:

Address 1:

Address 2 (optional):

City:

State/Province:

Postal code:

Country:

Phone:

Fax (optional):

Bill To

Same as order by:

Ship To

Same as order by:  Same as bill to:

Password:

Password confirmation:

Password must be at least 5 characters, and is case sensitive.

I accept the [Terms of Service](#):

Once you have completed entering in the required information, click on the **Register** button. You will be redirected back to the **Welcome** page.

## Confirmation

Upon completion of registration, the message below will appear.

Follow the instructions to verify your registration.

Upon verification, you will be directed to the Atlantis™ WebOrder **Welcome** page.

### Atlantis™ WebOrder

Your registration information has been received and a verification has been sent to the email you provided ([erik.carlsson@astratech.com](mailto:erik.carlsson@astratech.com)). Please look for this email and click on the verification link to confirm your registration. Reminder that all orders placed will be processed after approval of your registration.

<https://atlantisorder.astratechdental.com/>

## Orders

New orders are entered from the **Orders** page. Click the **Enter new order >>** link to begin.

You will be guided through the order process based on your specific selection.

Click the **Save and exit** button at any time to save your entries and return later to complete and submit your order to Astra Tech.

Go directly to any order by clicking on the hyperlink in the **Patient or PO/Ref.** column.

Sort any column by clicking on the headers.

Check shipping status of an order by clicking on the hyperlink in the **Status** column.

**Due Date** indicates date of delivery.

The screenshot shows the ATLANTIS WebOrder interface. On the left is a navigation menu with options: My Account, Preferences, Design Profiles, **Orders** (highlighted with a red box), Help, and Logout. The main content area is titled 'Orders' and features a search bar and an 'Advanced >' link. Below this is a table with the following data:

Patient or PO/Ref.	Prescriber	Units	Submit Date	Due Date	Status	Case Number	Delete
123345		4			Partial		X
Tally (this page): 1		4					

Below the table, the 'Enter new order >>' link is circled in red. A note states: 'All orders matching current query and with a Submit Date within selected range will be downloaded.' At the bottom, there is a date range selector showing '2009-03-01' to '2009-03-19' and a 'Download Data' button. The export file format is noted as CSV.

Entering an order

**Order page 1 – Overview** is used to enter initial case information.

Fields with an “\*” and implant location(s) are required.

**Process case – do not send images.** Select this option for your case to be fully processed. Abutment design images will not be sent.

**Hold case – send images.** Select this option to receive images of your case for review and approval prior to fabrication.

Indicate the **Model Type** being sent by clicking appropriate button.

Click the tooth number(s) to indicate where the implant(s) is located.

**Note:** The type of model you send for your Atlantis™ case will define the design options available to you.

*If sending a stone model of soft tissue, only an emergence width option of No Tissue Displacement will be available.*

*If sending a stone model with no soft tissue captured, "Use Subgingival Depth" margin position will not be available. The only emergence width option will be Full Anatomical Dimensions.*

Click **Next** to continue or **Save and exit** to save your progress and return later to complete your order.

**Order page 2 – Items** is used to specify implant and abutment information per tooth.

Select **Implant type** and size for each tooth number from the drop down menus. Select abutment material for each tooth number by clicking the corresponding button.

**Note:** Indicate if you want an Atlantis Gemini™ abutment by clicking the button in addition to the material selection.

**ATLANTIS™**  
WebOrder

English

Overview **Items** Parallelism Bill & Ship

My Account  
 Preferences  
 Design Profiles  
**Orders**  
 Help  
 Logout

Case Data & Design Preferences

**Important:** Please indicate specific abutment design preferences using the "edit >>" link(s) below. Otherwise, your abutment(s) will be designed using default design preferences.

**Note:** Provisional restorations, removable diagnostic preparations and positioning templates are not available for any orders containing Atlantis™ abutment in zirconia.

Tooth	Implant Type	Gold-Titanium			Atlantis Gemini™	Design Preferences
		shaded titanium	Zirconia	Zirconia		
		all	all	all	all	show all
45	Astra Tech	<input type="radio"/>	<input type="radio"/>	<input checked="" type="radio"/>	<input type="checkbox"/>	show edit>>
24	BioHorizons	<input type="radio"/>	<input checked="" type="radio"/>	<input type="radio"/>	<input type="checkbox"/>	show edit>>
21	Straumann	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="checkbox"/>	show edit>>
11	OsseoSpeed 4.5/5.0 (Lilac)	<input type="radio"/>	<input checked="" type="radio"/>	<input type="radio"/>	<input type="checkbox"/>	show edit>>
12	OsseoSpeed 4.5/5.0 (Lilac)	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="checkbox"/>	show edit>>

Additional Case Instructions  
 Please enter notes below in English.

Save and Exit Back Next


Click **Next** to continue.

Click **Save and exit** to save your progress and return later to complete your order.


Entering an order

To view **Design preferences** click the **show** link.

**Note: Healing abutment diameter** information is also entered on this page. Click the **hide** link to close the preferences.

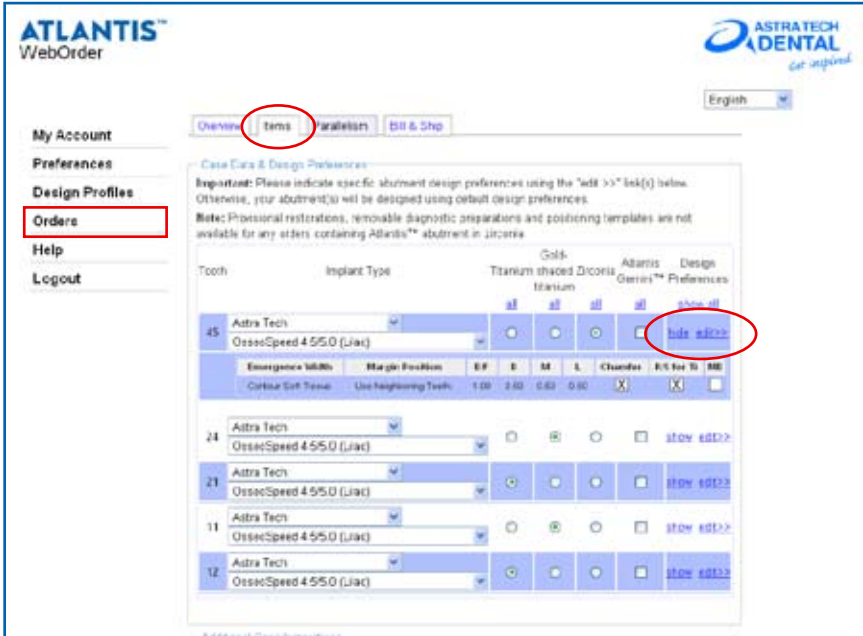

Click the **edit** link to change any preferences for each tooth number. Changes made here will not alter your default preferences. After making changes on the **Edit design preferences** page, click either the **Save changes** or **Cancel changes** button to exit.

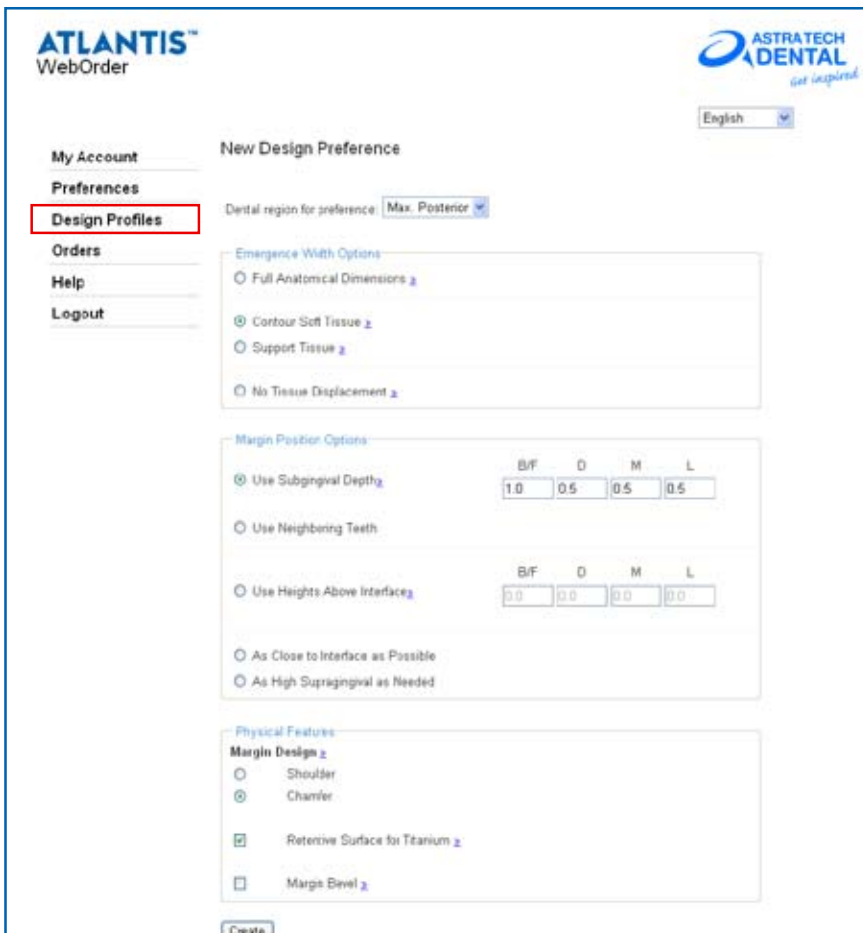

**Order page 3 – Parallelism** is used to indicate which abutments should be made parallel. Parallel abutments are indicated by the same color. Abutments that are to be parallel with each other should be checked in the same column and should be the same color on the chart (see example below).

**ATLANTIS™**  
WebOrder

ASTRATECH DENTAL  
Get inspired

English

Overview Items **Parallelism** Bill & Ship

My Account  
 Preferences  
 Design Profiles  
**Orders**  
 Help  
 Logout

Parallel Abutments

18 17 16 15 14 13 **12** **11** | **21** 22 23 **24** 25 26 27 28

48 47 46 **45** 44 43 42 41 | 31 32 33 34 35 36 37 38

If *final* restorations will be splinted the abutments MUST be designed to be parallel. Below, please select groups of abutments that should be designed parallel to one another. Parallel units should have the same "P" column selected and will have same background color in the arch diagram.

Abutment	P1	P2	P3	P4
12	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
11	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
21	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
24	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>

Abutment	P5
45	<input checked="" type="radio"/>

In the example below, abutments 11 and 21 are parallel. Abutment 12 is not parallel to 11, 21 or 22. Abutment 22 is not parallel to 12, 11, or 21.

Abutment	P1	P2	P3	P4
12	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>
11	<input type="radio"/>	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>
21	<input type="radio"/>	<input checked="" type="radio"/>	<input type="radio"/>	<input type="radio"/>
22	<input type="radio"/>	<input type="radio"/>	<input checked="" type="radio"/>	<input type="radio"/>

Save and Exit Back Next

© 2006-2009 Astra Tech. Patent #5674069, 6231342, 5989029. Additional patents pending. All trademarks are the property of their respective owners.

Click **Next** to continue.

Click **Save and exit** to save your progress and return later to complete your order.

**Order page 4 – Bill & ship** shows the remaining information populated from your previously entered preferences. To ensure accuracy, review all information in detail before submitting your order.

**ATLANTIS™**  
WebOrder

ASTRATECH  
DENTAL  
*Get inspired*

English

Overview Items Parallellism **Bill & Ship**

**My Account**

**Preferences**

**Design Profiles**

**Orders**

**Help**

**Logout**

Please review all shipping and billing information for accuracy before submitting your order.

**Payment method**  
Invoice

**Shipping preference**

- FedEx next business day, standard
- FedEx next business day, AM
- FedEx 2 day

Bill to this FedEx shipping account (XXXXXXXX):

**Bill To\***

Select address from list, or specify new: AstraTechDental

Business name: AstraTechDental

Address 1: Aminogatan 1

Address 2 (optional):

Address 3 (optional):

City: Mölndal

Postal code: 43121

Country: Sweden

Phone (optional): +46 317763012

Fax (optional):

**Ship To\***

Select address from list, or specify new: AstraTechDental

Business name: AstraTechDental

Address 1: Aminogatan 1

Address 2 (optional):

Address 3 (optional):

City: Mölndal

Postal code: 43121

Country: Sweden

Phone (optional): +46 317763012

Fax (optional):

Back Submit

© 2006-2009 Astra Tech. Patent #5674069, 6231342, 5989029. Additional patents pending. All trademarks are the property of their respective owners.

Click the **Submit** button to complete and send your order electronically to Astra Tech.

## Order ticket


Once your Atlantis™ WebOrder is submitted, an **Order ticket** containing a summary of your order information will appear. **Your Order ticket** identifies your case when it is received at Astra Tech.


Print your **Order ticket** by clicking the **Print this** link located at the bottom of the page. Include all items listed under **Case submission check list**.


This **Order ticket** must be signed, packaged with case materials and sent to Astra Tech in order for your case to be processed.

Click **View complete order information** to review all entries for this case.

Click the **Order list** link to return to the **Orders** page where you can check the status of any order entered through Atlantis™ WebOrder.





  
AWO0000065961

Heliungatan 1, SE-431 53 Mölndal, Sweden  
**Telefon +46 31 356 85 00**  
Email: [atlantis.europe@astratech.com](mailto:atlantis.europe@astratech.com)

**Patient or PO/Ref.:** 890+

**Prescriber:**  
**License:**  
**Phone:**

Tooth	Implant Type	Abutment Material	Atlantis Gemini™	Parallel Abutments
11	Astra Tech OsseoSpeed 3.0 (Yellow)	Gold-shaded titanium	<input type="checkbox"/>	P1

**Prescribing Clinician or Laboratory Signature:** *The signature below confirms that this product is being ordered at the request of a licensed dentist or on behalf of a licensed dentist whose information is on file with the laboratory.*

Please sign here

AstraTech AB/Marknad Dental  
Aminogatan 1  
Mölndal, 43121  
Sweden  
Phone: +46 317763012  
Fax:  
[christine.milkkrantz@astratech.com](mailto:christine.milkkrantz@astratech.com)

**Case submission check list**

To ensure efficiency in the processing of your order, the following materials are required:

- This Prescription Form
- Articulated Master Casts with analogs placed
- Bite Registration (recommended)
- Model of Soft Tissue

**NOTE:** Case material(s) must be disinfected prior to shipment.



[Print this](#) | [View complete order information](#) | [Order list](#)

[www.astratechdental.com](http://www.astratechdental.com)



## My account

On the **My account** page, the account manager can edit account information by clicking the **Edit this information >>** link and add new users by clicking the **Add new user >>** link.



English

---

**My Account**

- Preferences
- Design Profiles
- Orders
- Help
- Logout

---

**My Account**

Name & Email	Ordering Address
Tony Test1 tonytest1@hotmail.com	Tony Test Hagvagen Kiruna, 43525 Sweden Phone: +46 7763417 Fax: tonytest1@hotmail.com

[Edit this information >>](#)

**Account Users**

User	Email	Status	Delete
Tony Test1	tonytest1@hotmail.com	Acct. Mgr	

[Add new user >>](#)


You, as Account Manager, may:

- Add new users to your Atlantis™ WebOrder account
- Resend the verification email to Pending users
- Edit a user's name and password
- Delete users from your Atlantis™ WebOrder account

If you have any problems, please contact [Customer Service](#) at +46 31 356 85 00.

© 2006-2009 Astra Tech. Patent #5674069, 6231342, 5989029. Additional patents pending. All trademarks are the property of their respective owners.

Edit account information as needed and then click the **Save** button to return to the **My account** page.

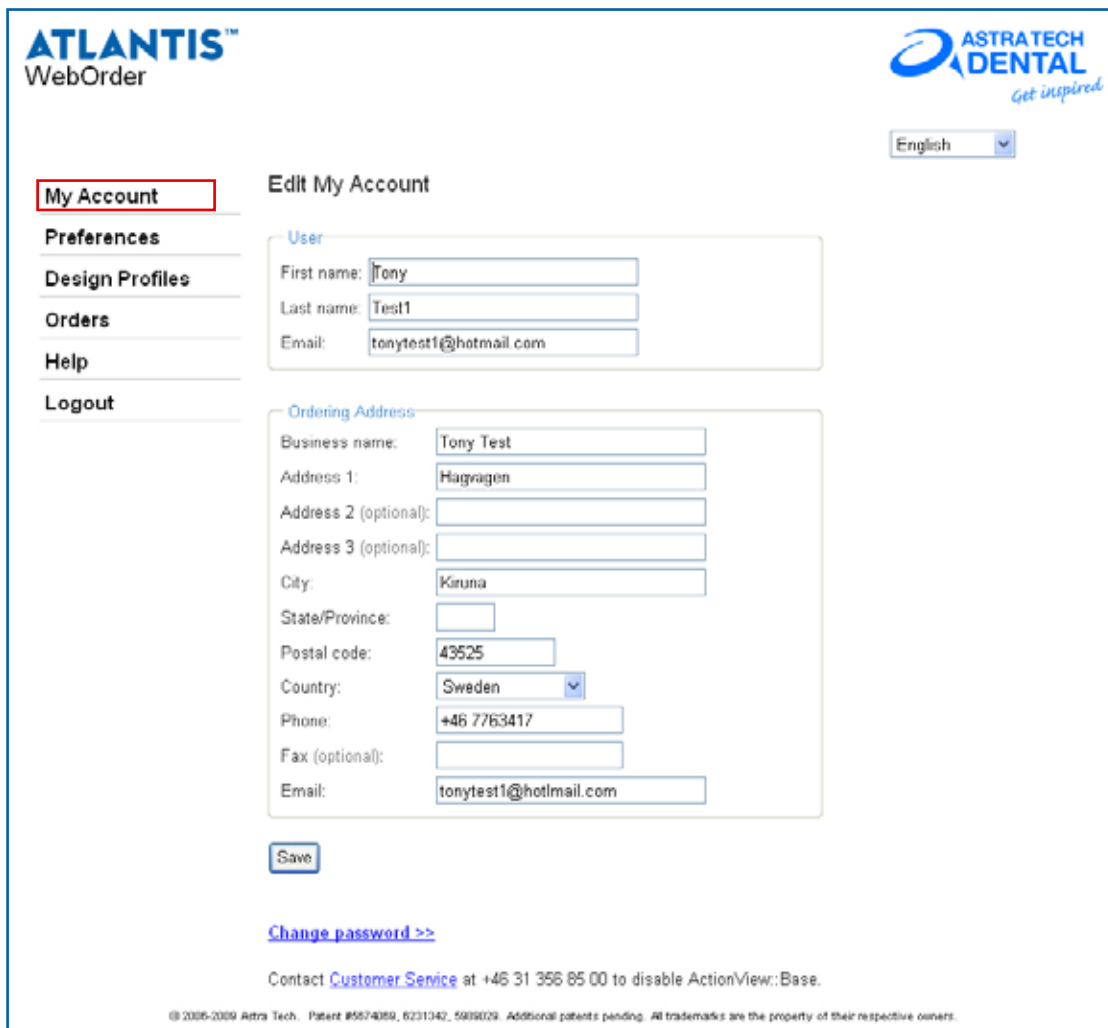

**ATLANTIS™**  
WebOrder

**ASTRATECH DENTAL**  
*Get inspired*

English

**My Account**

**Edit My Account**

**User**

First name:

Last name:

Email:

**Ordering Address**

Business name:

Address 1:

Address 2 (optional):

Address 3 (optional):

City:

State/Province:

Postal code:

Country:

Phone:

Fax (optional):

Email:

[Change password >>](#)


Contact [Customer Service](#) at +46 31 356 85 00 to disable ActionView:Base.

© 2005-2009 Astra Tech. Patent #5674089, 6231042, 5909029. Additional patents pending. All trademarks are the property of their respective owners.

To add a new user, enter the first name, last name and email address of the user. Click **Register** to add the new user. You will be redirected back to **My account** page.

The new user will receive an email to confirm their account.

Each user can access the same account but with different passwords, providing added flexibility in working with Atlantis™ WebOrder.

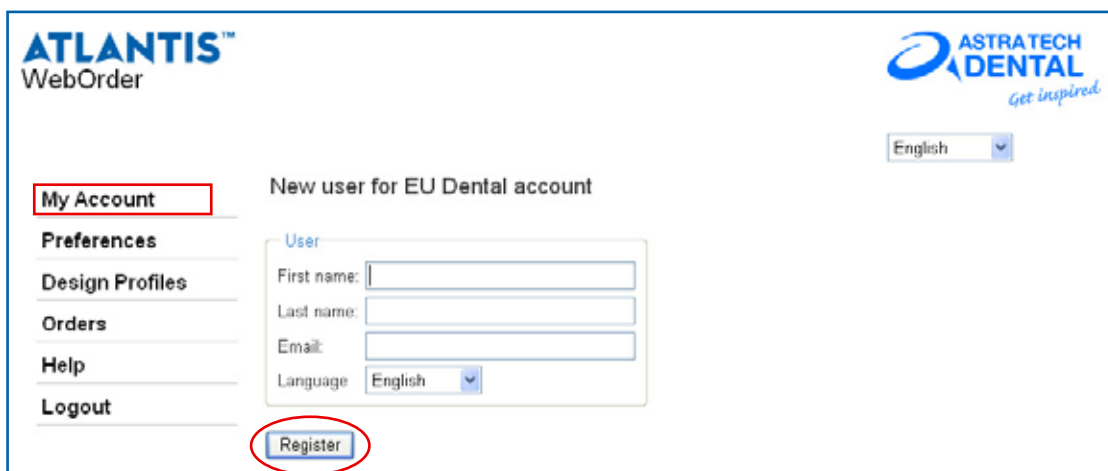

**ATLANTIS™**  
WebOrder

**ASTRATECH DENTAL**  
*Get inspired*

English

**My Account**

**New user for EU Dental account**

**User**

First name:

Last name:

Email:

Language:

# Preferences

Order preferences may be changed by clicking on the **Edit this information** link on the **Order Preferences** page.

**ATLANTIS™**  
WebOrder

**ASTRATECH DENTAL**  
Get Inspired

English

**My Account**  
[Preferences](#)  
[Design Profiles](#)  
[Orders](#)  
[Help](#)  
[Logout](#)

**EU Dental Account**  
**Order Preferences**

- Images - Hold case - send images. (eu.customer@intl.com)
- Abutment material - none
- Shipping carrier - Astra Tech selects the best way
- Shipping preference - end of business day
- Payment method - Invoice

**Saved Billing Addresses**

Business name	Address	Numbers	
EU Dental	Ansgatan 1 Mölndal, SE-431 21 Sweden	Phone: +46-31-7763000 Fac:	Default

**Saved Shipping Addresses**

Business name	Address	Numbers	
EU Dental	Ansgatan 1 Mölndal, SE-431 21 Sweden	Phone: +46-31-7763000 Fac:	Default

[Edit this information >>](#)

© 2007-2009 Astra Tech. Patent 850469, 821142, 848628. Additional patents pending. All trademarks are the property of their respective owners.

**ATLANTIS™**  
WebOrder

**ASTRATECH DENTAL**  
Get Inspired

English

**My Account**  
[Preferences](#)  
[Design Profiles](#)  
[Orders](#)  
[Help](#)  
[Logout](#)

**Edit Preferences**

**Images**

Process case - do not send images.  
 Hold case - send images.  
 Select per case.  
 Images sent to:

**Abutment material**

Titanium  
 Gold-shaded titanium  
 Zirconia  
 Select per case.

**Shipping carrier**  
 Astra Tech selects the best way  
 Please contact [Customer Service](#) at +46 31 356 86 00 to change this selection.

**Shipping preference**  
 Astra Tech selects the best way end of business day

**Payment method**  
 Invoice

**Billing Addresses**

Select address from list, or specify new: Tary Test

Business name:

Address 1:

Address 2 (optional):

Address 3 (optional):

City:

Postal code:

Country:

Phone (optional):

Fax (optional):

Set Default

**Shipping addresses**

Select address from list, or specify new: Tary Test

Business name:

Address 1:

Address 2 (optional):

Address 3 (optional):

City:

Postal code:

Country:

Phone (optional):

Fax (optional):

Set Default

**Save**


Once edited click **Save** to return to **Order Preferences**.


# Design profiles

In the **Design profiles** page you may add preferences by clicking on **add preference**.

You can also enter a new profile by clicking on **Enter new profile**.

Click **delete profile** to delete a specific profile.





English ▼

**My Account**

---

Preferences

---

Design Profiles

---

Orders

---

Help

---

Logout

### Design Profiles

**Default**

Dental Region	Emergence Width	Margin Position	B/F	D	M	L	Margin Design	Retentive Surface	Margin Bevel	Delete?
<a href="#">Default</a>	Contour Soft Tissue	Use Subgingival Depth	1.0	0.5	0.5	0.5	Chamfer	True	False	
<a href="#">Max. Anterior</a>	Full Anatomical Dimensions	Use Subgingival Depth	1.0	0.5	0.5	0.5	Chamfer	False	False	<a href="#">x</a>
<a href="#">Mand. Anterior</a>	Contour Soft Tissue	Use Subgingival Depth	1.0	0.5	0.5	0.5	Chamfer	True	False	<a href="#">x</a>
<a href="#">Mand. Posterior</a>	Contour Soft Tissue	Use Subgingival Depth	1.0	0.5	0.5	0.5	Chamfer	True	False	<a href="#">x</a>
<a href="#">add preference</a>										

**#007**

Dental Region	Emergence Width	Margin Position	B/F	D	M	L	Margin Design	Retentive Surface	Margin Bevel	Delete?
<a href="#">Default</a>	Full Anatomical Dimensions	Use Subgingival Depth	0.5	0.5	0.5	0.5	Chamfer	False	False	
<a href="#">add preference</a>										

[delete profile](#)

[Enter new profile >>](#)

© 2006-2009 Astra Tech. Patent #5674069, 6231342, 5989029. Additional patents pending. All trademarks are the property of their respective owners.

## Add preferences

On the **New Design Preference** page, you may edit and save your default for preferences such as **Emergence width options** (more information on page 22) and **Margin position options**. These preferences can be modified on an order-by-order basis.

The screenshot shows the 'New Design Preference' page in the ATLANTIS™ WebOrder system. The page is titled 'New Design Preference' and includes a sidebar with navigation options: My Account, Preferences, Design Profiles (highlighted), Orders, Help, and Logout. The main content area is divided into three sections:

- Dental region for preference:** A dropdown menu with options: Max. Anterior (selected), Max. Anterior, Mand. Anterior, and Max. Posterior.
- Emergence Width Options:** A section with a red circle around the title. It contains three radio button options:
  - Full Anatomical Dimensions
  - Contour Soft Tissue
  - Support Tissue
  - No Tissue Displacement
- Margin Position Options:** A section with two radio button options:
  - Use Subgingival Depth: Includes input fields for B/F (1.0), D (0.5), M (0.5), and L (0.5).
  - Use Neighboring Teeth
  - Use Heights Above Interface: Includes input fields for B/F (0.0), D (0.0), M (0.0), and L (0.0).
  - As Close to Interface as Possible
  - As High Supragingival as Needed
- Physical Features:** A section with a 'Margin Design' sub-section containing:
  - Shoulder
  - Chamfer
  - Retentive Surface for Titanium
  - Margin Bevel

A 'Create' button is located at the bottom of the form.

### Emergence width options

#### Full anatomical dimensions

Largest diameter abutment provided with the best emergence profile possible. A surgical incision for placement may be required.

#### Contour soft tissue

Medium diameter anatomically shaped abutment up to 1.0 mm larger than sulcus of model of the soft tissue model provided.

#### Support tissue

Anatomically shaped abutment that may be up to 0.2 mm larger than sulcus of the soft tissue model provided with desired emergence profile. "Easy" insertion.

#### No tissue displacement

Smallest diameter, non-anatomically shaped abutment with no soft tissue support. The abutment will not touch the soft tissue.

## Enter a new profile

A new design profile can be created using a customer name or other identifier. Selecting a specific design profile on the Overview page automatically uses the preferences specified for that profile. Click the link for “Enter new profile”, type in the profile name and click **Create**.





The screenshot shows the 'New Design Profile' page in the ATLANTIS WebOrder system. On the left, a navigation menu includes 'My Account', 'Preferences', 'Design Profiles' (highlighted with a red box), 'Orders', 'Help', and 'Logout'. The 'Design Profiles' section contains a 'Create' button circled in red. The main content area has a language dropdown set to 'English', followed by the title 'New Design Profile'. Below the title are two input fields: 'Design profile for prescriber:' and 'Prescriber postal code (optional):'.

After creating a new profile you may customize the profile by selecting the “Add preference” link under the Profile to access the **New Design Preference** page. (See previous page for details.)

The screenshot shows the 'New Design Preference' page. The left navigation menu is the same as in the previous screenshot, with 'Design Profiles' highlighted. The main content area has a language dropdown set to 'English' and the title 'New Design Preference'. Below the title is a dropdown menu for 'Dental region for preference:' with options: 'Max. Anterior', 'Max. Anterior', 'Mand. Anterior', and 'Max. Posterior'. The 'Emergence Width Options' section has radio buttons for 'Full Anatomical Dimensions', 'Contour Soft Tissue', 'Support Tissue', and 'No Tissue Displacement'. The 'Margin Position Options' section has radio buttons for 'Use Subgingival Depth', 'Use Neighboring Teeth', and 'Use Heights Above Interface', each with input fields for B/F, D, M, and L. The 'Physical Features' section has radio buttons for 'Margin Design' (Shoulder, Chamfer) and checkboxes for 'Retentive Surface for Titanium' and 'Margin Bevel'.

# Help functions

Visit the **Help** section for additional information related to ordering an Atlantis™ abutment.

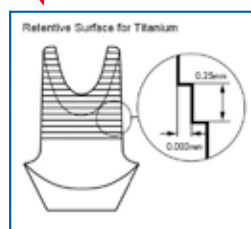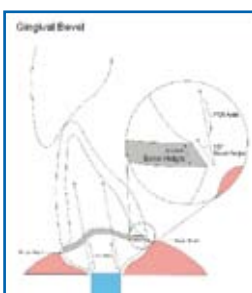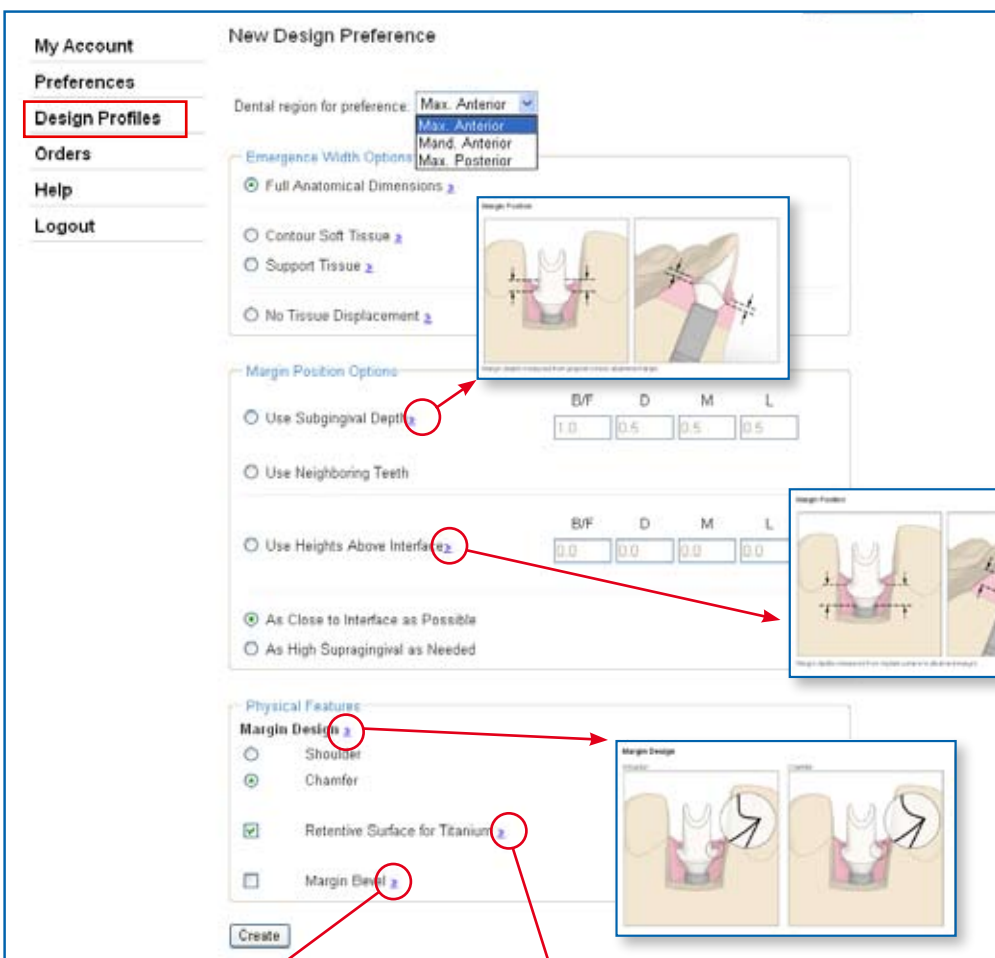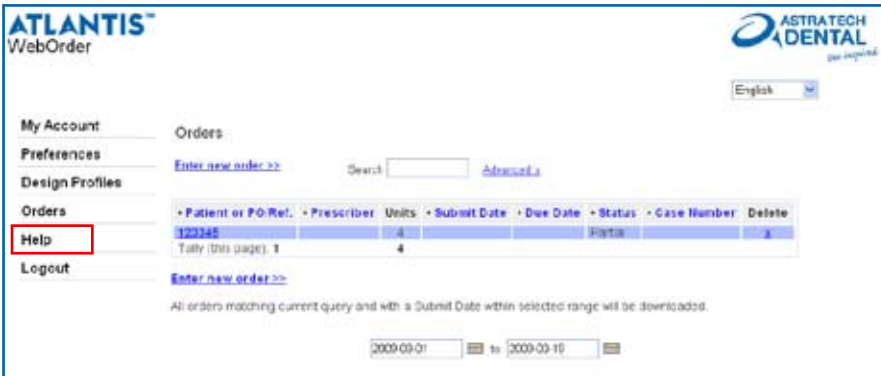

Clicking on any of the >> brings up an graphic illustration of the option.

## Emergence Width Options

**ATLANTIS™ WebOrder**

**My Account**  
 Preferences  
 Design Profiles  
 Orders  
 Help  
 Logout

**New Design Preference**

Dental region for preference: **Max. Posterior**

**Emergence Width Options**

- Full Anatomical Dimensions >>
- Contour Soft Tissue >>
- Support Tissue >>
- No Tissue Displacement >>

**Margin Position Options**

Use Subgingival Depth

BIF	D	M	L
1.0	0.5	0.6	0.6

Use Neighboring Teeth

Use Heights Above Interface

BIF	D	M	L
0.0	0.0	0.0	0.0

As Close to Interface as Possible

As High Supragingival as Needed

**Physical Features**

**Margin Design**

- Shoulder
- Chamfer

Retentive Surface for Titanium >>

Margin Bevel >>

**Create**

**Full Anatomical Dimensions**  
 Largest amount of abutment width and emergence profile possible in largest section for given body height for abutment.

**Contour Soft Tissue**  
 Reduces emergence and abutment width up to 1.0 mm above the neck of abutment without position.

**Support Tissue**  
 Reduces abutment width on up to 0.5 mm above the crown neck with tissue height provided and case given profile. Tissue position.

**No Tissue Displacement**  
 Reduces emergence, subgingival depth and neck with no soft tissue support. The abutment will be left alone.

Clicking on any of the >> brings up an graphic illustration of the option. See "Design guide – Atlantis™ patient-specific abutments" manual for more information.

## Model Type

**ATLANTIS™ WebOrder**

**ASTRATECH ADENTAL**  
*Get Inspired*

English

**My Account**  
 Overview | Items | Parallelism | B&I & Slip

**Preferences**

**Design Profiles**

**Orders**

**Help**

**Logout**

**Model Type**

- Process case - do not send images.
- Hold case - send images. (Case held until approval of abutment design images is provided)

Email for images (request multiple email addresses by a comma, 60 characters max.)  
 tonytaft@hotmail.com

**Model Type**

- Stose model with removable soft tissue >>
- Stose model of soft tissue >>
- Stose model - no soft tissue >>

**Case Overview**  
 Please indicate the implant locations.

18	17	16	15	14	13	12	11	10	9	8	7	6	5	4	3	2	1
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
46	47	45	44	43	42	41	31	32	33	34	35	36	37	38	39	40	48
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>

**Reminder:** For cases with four or more adjacent edentulous spaces a removable diagnostic preparation that relates to the master cast is required.

**Save and Exit** **Next**

**Removable Soft Tissue**  
 Reduces emergence and abutment width up to 1.0 mm above the neck of abutment without position.

**Support Tissue**  
 Reduces abutment width on up to 0.5 mm above the crown neck with tissue height provided and case given profile. Tissue position.

**No Tissue Displacement**  
 Reduces emergence, subgingival depth and neck with no soft tissue support. The abutment will be left alone.



# Atlantis™ – patient-specific abutments for all major implant systems



## Australia

Astra Tech Pty Ltd.  
Suite 1, 53 Grandview St, Pymble NSW 2073  
☎ +61 2 9488 3500. 📠 +61 2 9440 0744  
[www.astratechdental.com.au](http://www.astratechdental.com.au)

## Austria

Astra Tech GesmbH  
Schloßhofer Straße 4/4/19, AT-1210 Wien  
☎ +43-(0)1-2146150. 📠 +43-(0)1-2146167  
[www.astratechdental.at](http://www.astratechdental.at)

## Benelux

Astra Tech Benelux B.V.  
Signaalrood 55, NL-2718 SG Zoetermeer  
☎ +31 79 360 1955/+32 3 232 81 50  
📠 +31 79 362 3748/ +32 3 213 30 66  
[www.astratechdental.nl](http://www.astratechdental.nl)

## Canada

Astra Tech Inc.  
2425 Matheson Blvd East, 8th Floor  
Mississauga, ON L4W 5K4  
☎ +1 905 361 2844  
[www.astratechdental.com](http://www.astratechdental.com)

## Denmark

Astra Tech A/S  
Roskildevej 163, 1. th., DK-2620 Albertslund  
☎ +45 43 71 33 77. 📠 +45 43 71 78 65  
[www.astratechdental.dk](http://www.astratechdental.dk)

## Finland

Astra Tech Oy  
PL 96, FI-02231 Espoo  
☎ +358 9 8676 1626. 📠 +358 9 804 4128  
[www.astratechdental.fi](http://www.astratechdental.fi)

## France

Astra Tech France  
7, rue Eugène et Armand Peugeot, TSA 90002  
FR-92563 Reuil Malmaison Cedex  
☎ +33 1 41 39 02 40. 📠 +33 1 41 39 02 44  
[www.astratech.fr](http://www.astratech.fr)

## Germany

Astra Tech GmbH  
An der kleinen Seite 8. DE-65604 Elz  
☎ +49 6431 9869 0. 📠 +49 6431 9869 500  
[www.astratechdental.de](http://www.astratechdental.de)

## Italy

Astra Tech S.p.A.  
Via Cristoni, 86, IT-40033 Casalecchio di Reno (BO)  
☎ +39 051 29 87 511. 📠 +39 051 29 87 580  
[www.astratechdental.it](http://www.astratechdental.it)

## Japan

Astra Tech K.K.  
1-7-16 Sendagaya, Shibuya-ku, Tokyo 151-0051  
☎ +81 3 5775 0515. 📠 +81 3 5775 0571  
[www.astratech.jp](http://www.astratech.jp)

## Norway

Astra Tech AS  
Postboks 160, NO-1471 Lørenskog  
☎ +47 67 92 05 50. 📠 +47 67 92 05 60  
[www.astratechdental.no](http://www.astratechdental.no)

## Poland

Astra Tech Sp.z o.o.  
ul. Orężna 58 02-937 Warszawa  
☎ +48 22 853 67 06. 📠 +48 22 853 67 10  
[www.astratechdental.pl](http://www.astratechdental.pl)

## Spain

Astra Tech S.A.  
Calle Ciencias nº 73 derecha. Nave 9,  
Polígono Industrial Pedrosa,  
ES-08908 L'Hospitalet de Llobregat  
☎ Servicio al cliente: +34.902.101.558  
☎ +34.932.643.560. 📠 +34.933.363.231  
[www.astratechdental.es](http://www.astratechdental.es)

## Sweden

Astra Tech AB  
P.O. Box 14, SE-431 21 Mölndal  
☎ +46 31 776 30 00. 📠 +46 31 776 30 17  
[www.astratechdental.se](http://www.astratechdental.se)

## Switzerland

Astra Tech SA  
Avenue de Sévelin 18, P.O. Box 54  
CH-1000 Lausanne 20  
☎ +41 21 620 02 30. 📠 +41 21 620 02 31  
[www.astratech.ch](http://www.astratech.ch)

## United Kingdom

Astra Tech Ltd.  
Brunel Way, Stonehouse, Glos GL10 3SX  
☎ +44 1453 791763. 📠 +44 1453 791001  
[www.astratechdental.co.uk](http://www.astratechdental.co.uk)

## USA

Astra Tech Inc.  
590 Lincoln Street, Waltham, MA 02451  
☎ +1-800-531-3481. 📠 +1-781-890-6808  
[www.astratechdental.com](http://www.astratechdental.com)

## Other Markets

Astra Tech AB, Export Department  
P.O. Box 14, SE-431 21 Mölndal, Sweden  
☎ +46 31 776 30 00. 📠 +46 31 776 30 23  
[www.astratechdental.com](http://www.astratechdental.com)

